

### Rotary evaporator, Rotavapor® R-3



For both teaching and general purpose rotary evaporation applications.

- Simple, manual lift to raise and lower the evaporating flask
- Continuously adjustable rotation speed between 20rpm and 280rpm for high distillation rates
- Digital bath temperature display, easy to set and read
- Heating bath made from durable stainless steel
- Heating bath can accommodate up to 4L evaporating flasks
- Three function Combi-Clip included:
  - Secures flask to the vapour duct
  - Removes the flask from the vapour duct
  - Removes the vapour duct from the R-3
- Only one power cable is required to operate the R-3

The R-3 evaporator is easy to assemble and operate.

#### Technical Specification - General

Rotational speed, rpm	20 to 280
Flask size, mL	50 to 4,000
Flask load capacity, kg	3 (maximum)
Temperature range, °C	20 to 95
Capacity, L	4
Mass, kg	15
Electrical supply	220-240V, 50/60Hz, 60W

Catalogue No	Legacy No	Description
10201395	-	Rotavapor® rotary evaporator R-3


#### Accessories

Catalogue No	Legacy No	Alt. No	Description
10474252	-	48052	Protective shield against splashing or possible flask breakage
12752315	-	000636	Vacuum seal
12787173	-	48618	Flask support ring
10128981	-	36576	Reitmeyer adapter
10493682	-	71110	V-700 vacuum pump




## Have you thought about.....

### Thermometers, yellow backed, mercury filled, total immersion, Fisherbrand®

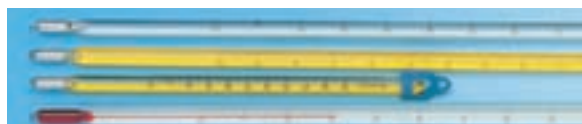

General purpose/stirring thermometers with amber graduations, reinforced bulb for stirring and non-roll fitting.

Catalogue No	Catalogue No	Range, °C	Subdivisions, °C	Length, mm
10080331	-	-10 to 50	0.5	155
10746761	-	-10 to 110	1	155
10684002	-	-10 to 250	5	155
10141851	-	-10 to 360	5	155

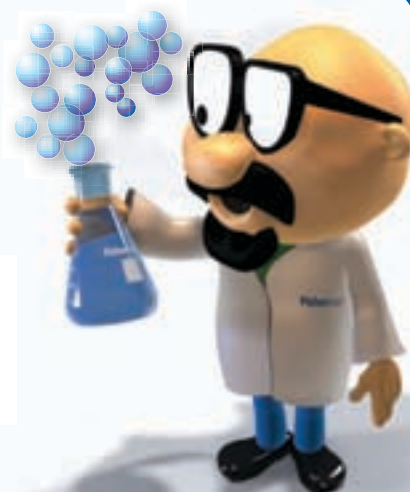

Look for the Fisherbrand® logo throughout the catalogue. It's your guide to great value products.