



Pipettors, variable volume, multichannel, Fisherbrand® Elite™



Fisherbrand®

- Fully autoclavable
- Extremely low plunger forces - work longer without fatigue and reduce the risk of RSI
- Soft touch tip ejection - a gearing mechanism reducing tip ejection force by up to 50%
- Ergonomic comfortable lightweight handle
- Optimal fit with Fisherbrand® SureOne™ pipettor tips
- Super blow out piston - creates 150% bigger air push than conventional piston ensuring delivery of micro size drops

8 channel

Catalogue No	Legacy No	Volume range, µL	Increments, µL	Accuracy, %	Precision, %
11815772	-	1 to 10	0.02	2.4 to 12.0	1.6 to 8.0
11835772	-	5 to 50	0.1	1.5 to 5.0	0.7 to 2.0
11825772	-	10 to 100	0.2	1.3 to 5.0	0.5 to 2.0
11845772	-	30 to 300	1	1.0 to 4.6	0.3 to 1.5

12 channel


Catalogue No	Legacy No	Volume range, µL	Increments, µL	Accuracy, %	Precision, %
11855772	-	1 to 10	0.02	2.4 to 12.0	1.6 to 8.0
11875772	-	5 to 50	0.1	1.5 to 5.0	0.7 to 2.0
11865772	-	10 to 100	0.2	1.3 to 5.0	0.5 to 2.0
11885772	-	30 to 300	1	1.0 to 4.6	0.3 to 1.5



Bulk pack



Racked open

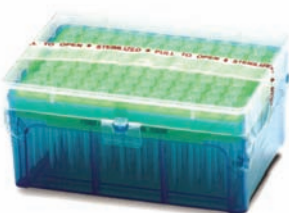

Racked sterile

Pipettor tips, standard, universal fit, Fisherbrand® SureOne™



Fisherbrand®

Fisherbrand® SureOne™ standard pipettor tips are a comprehensive line of universal fit pipettor tips offered in many packaging configurations, including the environmentally friendly reload. Tips are available in a range of volumes, from 5µL to 10mL, and are compatible with most popular brands of pipettor (optimal fit with Fisherbrand® Elite™ pipettor).

SureOne™ insert

- Contact clip prevents unnecessary contact with the tips while clipping in inserts
- Lot numbered inserts for traceability include tip volume
- Indented insert holes inhibit tip movement, stabilising them for multichannel work
- Rigid inserts are all 96 well multichannel format, offering a robust non-pliable design providing secure seals on every tip

SureOne™ rack

- Flush fit inserts sit flush with rack, allowing easy multichannel pipetting
- Recycle logo moulded into all rack components, allows recycling companies to properly recycle the plastic
- Notched clasp designed to prop lid open for autoclave venting
- SureStack™ lid rim prevents racks from sliding around when stacked and moved around the laboratory
- Safe lock latch allows for one handed opening and closing ability and will not open if dropped
- Sure Snap design securely holds insert in place, preventing it from being lifted out while pipetting
- Flip lid can be either a hinged or lift off lid (hinged format does not prevent use of a multichannel pipettor)

SureOne™ reload packaging

- Product information displayed on box
- Dust flap prevents debris from contaminating the product
- Finger knockouts allow easy access to slide out product
- Canyon cut package allows for easy tip insert removal