

### Pipettors, electronic, Picus



Ideal for use in pipetting, reverse pipetting, manual pipetting, multi-pipetting, dilution, sequential dispensing, multi-aspiration and titrate applications.

- Ergonomic lightweight design to protect user from Repetitive Strain Injury (RSI)
- Electronic brake and piston control system guarantee accurate and precise pipetting
- Intuitive user interface
- Fast volume setting and menu navigation
- Hot key for quick access to favourite programs
- Calibration in 1, 2 or 3 points
- Autoclavable lower parts
- Li-Polymer battery allows charging in approx. 1 hour

#### Technical Specification - Specific


##### Single channel

	13456358	13466358	13486358	13476358	13496358	13406368
Volume range, $\mu\text{L}$	0.2 to 10	5 to 120	10 to 300	50 to 1,000	100 to 5,000	500 to 10,000
Volume, test, $\mu\text{L}$	1 to 10	12 to 120	30 to 300	100 to 1,000	500 to 5,000	1,000 to 10,000
Inaccuracy, %	0.90 to 2.50	0.40 to 2.00	0.40 to 1.50	0.40 to 1.50	0.50 to 1.00	0.60 to 3.00
Imprecision, %	0.40 to 1.50	0.40 to 2.00	0.15 to 1.00	0.15 to 0.50	0.15 to 0.40	0.20 to 0.60

##### Multichannel

	13416368	13426368	13436368	13446368
Volume range, $\mu\text{L}$	0.2 to 10	5 to 120	10 to 300	50 to 1,200
Volume, test, $\mu\text{L}$	1 to 10	12 to 120	30 to 300	120 to 1,200
Inaccuracy, %	0.90 to 4.00	0.50 to 2.00	0.50 to 2.00	0.50 to 2.50
Imprecision, %	0.50 to 3.00	0.20 to 1.50	0.20 to 1.00	0.20 to 1.00

	13456368	13466368	13476368	13486368
Volume range, $\mu\text{L}$	0.2 to 10	5 to 120	10 to 300	50 to 1,200
Volume, test, $\mu\text{L}$	1 to 10	12 to 120	30 to 300	120 to 1,200
Inaccuracy, %	0.90 to 4.00	0.50 to 2.00	0.50 to 2.00	0.50 to 2.50
Imprecision, %	0.50 to 3.00	0.20 to 1.50	0.20 to 1.00	0.20 to 1.00

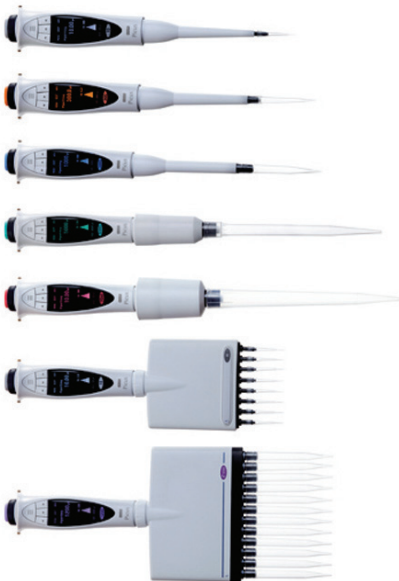

##### Single channel

Catalogue No	Legacy No	Alt. No	Volume range, $\mu\text{L}$	Increments, $\mu\text{L}$
13456358	-	735021	0.2 to 10	0.02
13466358	-	735041	5 to 120	0.10
13486358	-	735061	10 to 300	0.20
13476358	-	735081	50 to 1,000	1.00
13496358	-	735101	100 to 5,000	5.00
13406368	-	735111	500 to 10,000	10.00

##### 8 channel

Catalogue No	Legacy No	Alt. No	Volume range, $\mu\text{L}$	Increments, $\mu\text{L}$
13416368	-	735321	0.2 to 10	0.02
13426368	-	735341	5 to 120	0.10
13436368	-	735361	10 to 300	0.20
13446368	-	735391	50 to 1,200	1.00

##### 12 channel

Catalogue No	Legacy No	Alt. No	Volume range, $\mu\text{L}$	Increments, $\mu\text{L}$
13456368	-	735421	0.2 to 10	0.02
13466368	-	735441	5 to 120	0.10
13476368	-	735461	10 to 300	0.20
13486368	-	735491	50 to 1,200	1.00

##### Accessories

Catalogue No	Legacy No	Alt. No	Description
13104254	-	730981	Charging stand for one pipettor
13157544	-	730991	Charging carousel for four pipettors
13194424	-	725620	Linear stand for all Sartorius Biohit pipettors