


Funnels, sample chamber, disposable, Thermo Scientific Shandon EZ Cytofunnel®

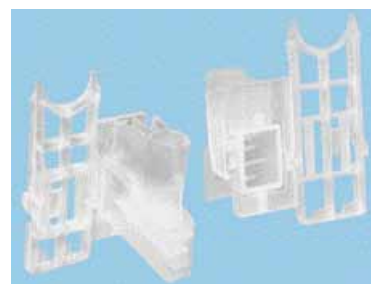

A range of cytofunnels intended for single use with all Thermo Scientific Shandon Cytospin® cytocentrifuge applications.

- Optimise cell yield
- Safe – reduce the risk of exposure to pathogenic samples
- Disposable – eliminate time-consuming decontamination and cleaning
- Decrease the possibility of cross-contamination
- Eliminate the need for metal Cytoclips®
- Closure caps provided with kit
- Produce more cost-effective thin-layer preparations when compared to other thin-layer methods
- Simple and easy to use operation illustrations are included in the package

EZ Megafunnel® with caps and slides

Designed for convenient preparation of larger volume samples – up to 6mL, minimising the need to pre-centrifuge. Can prepare both air-dried and fixed preparations, providing additional valuable diagnostic information. Megafunnel Cytoslides® included.

Cell deposition area of 22mm x 14.75mm (325mm²) for increased cell yield.

EZ Single Cytofunnel® with white filter cards

Use for sample volumes up to 0.5mL. Pre-attached white filter card for high absorbency.

Cell deposition area of 6mm in diameter (28mm²).

EZ Single Cytofunnel® with brown filter cards

Use for sample volumes of 0.4mL or less. Pre-attached brown filter card for slower absorbency. Excellent for scanty specimens such as CSF.

Cell deposition area of 6mm in diameter (28mm²).

EZ Double Cytofunnel® with white filter cards and caps

Allows two sample deposition areas on a single slide: double the number of sample presentations per run. Funnel markings correspond to slide markings for simple identification. Pre-attached filter card.

Catalogue No	Alt. No	Description	Pack qty
SDA78710001	A78710001	EZ Megafunnel® with caps and slides	25
SDA78710003	A78710003	EZ Single Cytofunnel® with white filter cards and caps	40
SDA78710020	A78710020	EZ Single Cytofunnel® with white filter cards and caps	200
SDA78710004	A78710004	EZ Single Cytofunnel® with brown filter cards and caps	40
SDA78710005	A78710005	EZ Double Cytofunnel® with white filter cards and caps	40
SDA78710021	A78710021	EZ Double Cytofunnel® with white filter cards and caps	200

Accessories

Catalogue No	Alt. No	Description	Pack qty
SD5991023	5991023	Brown filter cards	200
SD5991039E	5991039	Double white Cytofunnel® with disposable sample chamber	25
SD5991040	5991040	Single white Cytofunnel® with white filter card and caps	25
SD5991022	5991022	TPX white filter cards	200

All products are for research use only, unless otherwise indicated.